



Field Service Procedure

Part Number: SP00290
Rev: -
Date: 25 September 2003
© 2003 Draeger Medical, Inc.

GCX Pivot Rail Installation

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

GCX Mounting Assembly Installation Guide

Purpose: The purpose of this document is to provide instruction for the mounting of the Pivot Rail Assembly to Dräger Medical Anesthesia Machines.

Parts Reference

The following parts and hardware are included with this installation kit.
(see illustrations for parts, hardware not shown):

Item #	Description	Qty
1	Pivot Rail Assembly	1

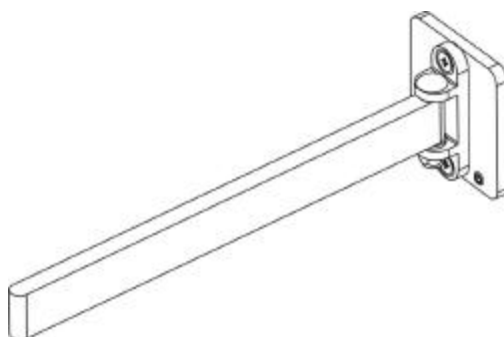


Figure 1. GCX Pivot Rail Assembly

Tools required for this procedure: 3mm hex key (provided).

INSTALLATION PROCEDURE:

Mounting the Pivot Rail Assembly (left side)

1. Insert the Slide (1) in the GCX channel (with clamping screws at bottom) and move the Pivot Rail Assembly to the desired height.
2. Using the 3mm hex key (provided), tighten the two (2) clamping screws to secure the assembly in the channel.

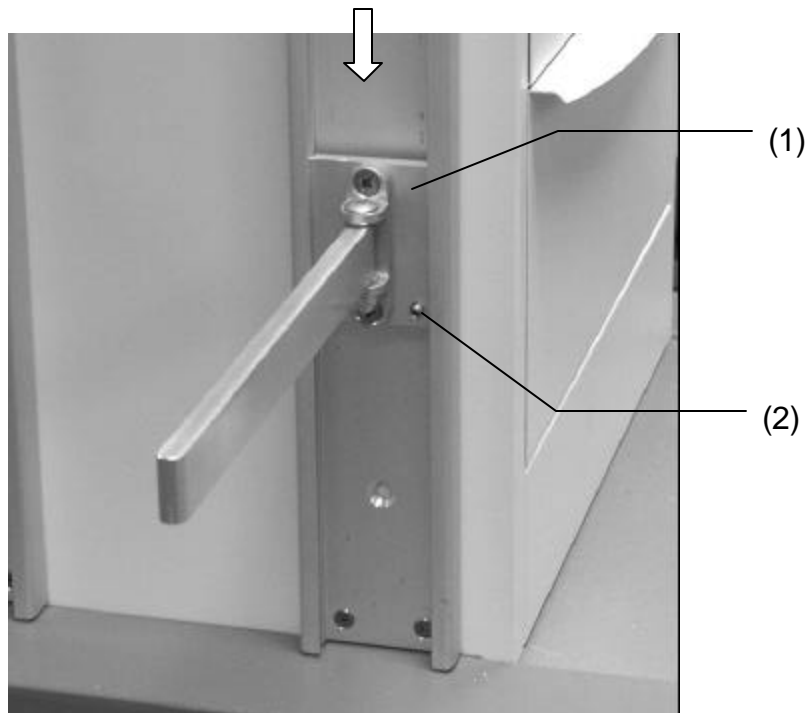


Figure 2. Pivot Rail Assembly

INSTALLATION PROCEDURE:

Mounting the Pivot Rail Assembly (right side)

1. Insert the Slide (1) in the GCX channel (with clamping screws at bottom) and move the Pivot Rail Assembly to the desired height.
2. Using the 3mm hex key (provided), tighten the two (2) clamping screws to secure the assembly in the channel.

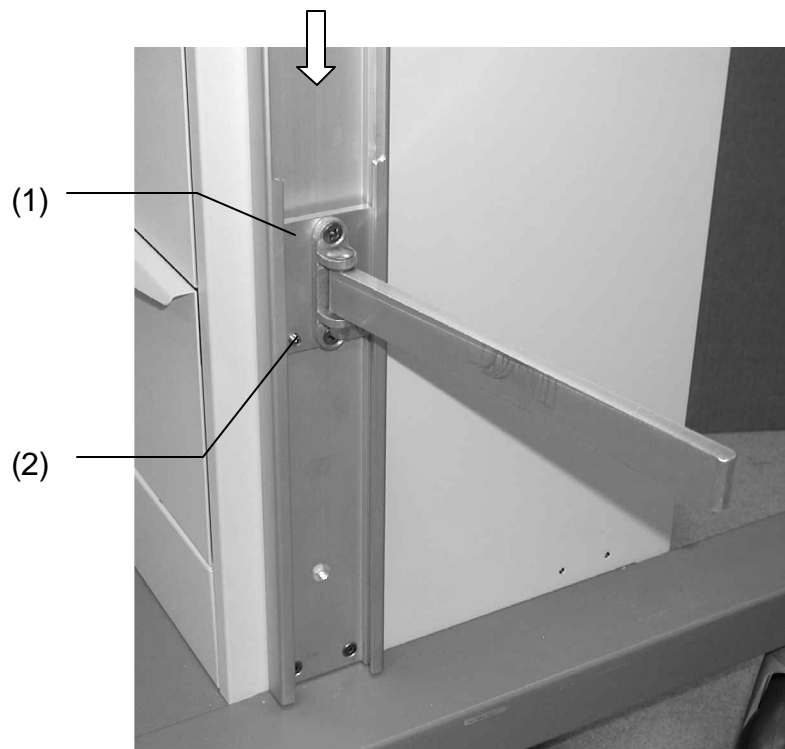


Figure 3. Pivot Rail Assembly

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)



DrägerService is a division of
Draeger Medical, Inc.
3122 Commerce Drive
Telford, PA 18969
Tel: (215) 721-5402
(800) 543-5047
Fax: (215) 721-5784
Web: www.draegermedical.com
Printed in the U.S.A.